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User Manual

BPB-M02H

Enjoy it.

V1.0
908

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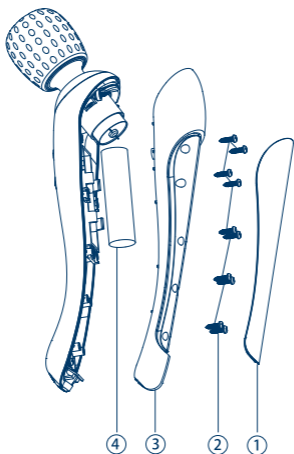
Features

1. Waterproof effect.
2. Five different vibration modes are available , which is chosen to suit individual comfort .
3. Use LED indicator, buzzer prompt, more convenient to use.
4. USB fast charging, built-in large-capacity lithium battery power supply, rechargeable, easy to use.
5. Intelligent vibration mode memory function and button lock function.
6. Built-in 15 minutes timed automatic shutdown function.

Instructions

1. Please fully charged before using this product for the first time.
2. When the battery power is insufficient, the LED indicator light will go into the flickering state with " ticking " prompt sound, the product will automatically shut down 10 seconds later.
3. During charging operation, connect the charging cable interface to the charging terminal of the product body , and then plug the USB port of the other terminal of the charging cable into the mobile power or computer USB port (5 V DC power supply) to charge, meanwhile, the charging indicator light is in the stage of flashing.

Operation instructions

1. Press the "⏻" button shortly to turn on the product, the product body will start to enter the vibration massage mode, the timing time is 15 minutes, and the work will stop automatically after 15 minutes.
2. Then press shortly the "🔄" button to select five vibration modes of "1 (weak), 2 (moderate), 3 (strong), 4 (intermittent mode), 5 (wave mode)" , the corresponding LED indicator light will blink"once , twice, three times, four times, five times. Cycle in order : 2-3-4-5-1-2-3-..... (After turning off the power and then turn it on , it will automatically return to the vibration , massage mode of the previous shutdown.

Instructions

3. Press shortly the “⏻” button again and the product body will enter the shutdown state.
4. After turning off the power, press the “⏻” button and “🔄” button simultaneously for about 2 seconds with two “ticking” prompt sound from the product to activate button lock function. If you want to unlock, press the “⏻” and “🔄” button simultaneously for about 2 seconds, deblocking achieves after “ticking” prompt sound occurs.


Use of rechargeable battery

1. During the use, if the charging indicator light is always flashing with “ticking” prompt sound, which means the built-in lithium battery is low, insert the DC port of the charging cable into the product body, and then plug the USB port at the other terminal of the charging cable into mobile power or computer USB port (5V DC power supply) to get charged.
2. When charging, the charging indicator light of product body is continuously flashing. After the battery is fully charged, the charging indicator light is always on. It takes about 3-4 hours to be charged fully. Work continuously about 1-2 hours after charged fully.
3. If you do not charge the battery for a long time, it will shorten the service life of the rechargeable battery. It is recommended to charge the battery fully at least once a month.

Attentions

1. Please read carefully before using this product, understand the full content of user manual before operation.
2. Prohibition use in pregnant women, children, heart disease, abnormal blood pressure, malignant tumors, cerebrovascular patients, acute patients or other patients accepting treatment of doctor.
3. Those who have no ability to express their consciousness are forbidden to use.
4. Patients with skin sensory disturbances are forbidden to use.
5. No use within 1 hour after meal.
6. If you feel unwell or have abnormal skin during use, please stop using it immediately and consult a doctor.
7. It is forbidden to use this product in high heat, flammable, electromagnetic radiation and wet places.
8. It is forbidden to disassemble this product privately to avoid malfunction resulting from collision and water ingress.
9. It is recommended that use time should not exceed 15 minutes (the product will automatically shut down after 15 minutes). If you need to continue using it, please start the product again after 30 minutes.
10. Excessive continuous use time may cause product overheating and shorten the service life. If the product is overheated, please stop using it and you can use it again until the product cools down.
11. When this product is discarded, please deal with it in accordance with the relevant environmental protection regulations, otherwise it may cause environmental pollution.

Specifications

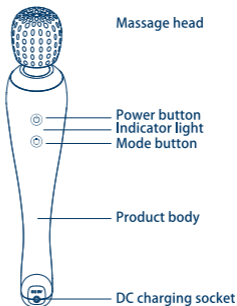Product specifications

Product size: 285×58mm

Net weight : about 360g

Power: DC 5V/1A (USB)

Input power: 5W

Lithium battery : DC 3.7V/2200mAh

Charging time: about 3-4 hours

Normal working conditions :

Ambient temperature: +5°C~+40°C

Relative humidity: ≤ 80%

Atmospheric pressure: 86kPa~106kPa

Transportation and storage conditions

Ambient temperature range: -20°C~+55°C

Relative humidity range: ≤ 93%

Atmospheric pressure range: 50kPa~106kPa

Cleaning, maintenance and storage

1. Cleaning and maintenance

- 1) Please clean the product with a neutral detergent. After cleaning, wipe it off with a soft cotton cloth.
- 2) Do not use any wire mesh, iron brush, sharp tools etc., and clean it with volatile, corrosive solvents.
- 3) Usually, when not in use, put the product in a box and store it in a dry and cool place.
- 4) Do not store this product near high temperatures or open flames, or in direct sunlight, which may cause color fading and damage.
- 5) Avoid using corrosive liquids such as benzene, thinner and other chemical solvents on this product and accessories.

2. Storage

- 1) Please put it in a place where children can't touch.
- 2) Do not keep it in direct sunlight, high temperature and wet places.
- 3) Please keep it in a dry and ventilated place.
- 4) Please do not disassemble, repair or modify this product privately, which may cause an accident.

Do not disassemble or try to repair privately.

Faults:

1. The product can not be started.
2. Can't charge or no charging display.

Inspection matters:

1. Is the internal rechargeable battery low? Whether the product is in the button lock mode.
2. Is the USB wire connected well ?
3. Is the charging power supply low?

Solutions:

1. If the internal rechargeable battery is low, please charge it with the USB cable. About button unlocking, please follow the instructions to unlock it.
2. Check if the plugs on both ends of the USB are in place
3. Check if the charging power supply is normal.

Troubleshooting


Environmental protection and resource recycling

Lithium batteries are used in this product. Make sure the old batteries are handled uniformly at designated locations in accordance with local regulations.

The way to remove the built-in rechargeable battery

This product contains a rechargeable battery, before the consumer decides to discard the product, be sure to remove the battery from the product.

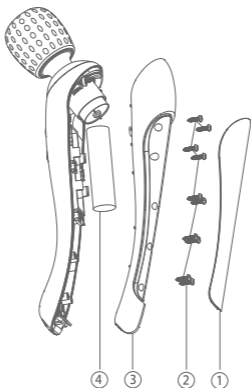1. Must turn off the power before you decompose the product and remove the battery.
2. Remove the battery according to the sequence steps 1 to 4 shown in the right diagram.
3. The scrapped battery must be disposed of in a government-designated place or recycling organization.

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